Hand-held massager

Massageapparat / Massasjeapparat /
Massagegerät / Hierontalaite

ENG

SE

NO

DE

FI

Thank you for choosing to purchase a product from Rusta!

Read through the entire manual before installation and use!

Hand-held massager

Read the user manual carefully to make sure that the product is assembled, used and maintained correctly, as described in the instructions in this user manual. Keep the user manual for future reference.

SAFETY INSTRUCTIONS

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be done by children unless they are older than 8 and supervised. Keep the appliance and its cord out of reach of children less than 8 years.
- Keep the massager away from excessively humid, dusty environment and water.
- Please do not use the product within one hour before or after meals or drinking.
- Please do not stress or stand on the massager to avoid damage.
- Please do not use it when drowsy or sleeping.
- Close supervision is necessary when this appliance is used by children.
- Do not use on the areas where skin is swollen, bruised, or inflamed.
- If the supply cord is damaged, it should be replaced by supplier, retailer, it's service agent or a similarly qualified person in order to avoid a hazard.


WARNING!

If these instructions are not followed, an accident may lead to serious personal injury or property damage.


PRODUCT

MASSAGER

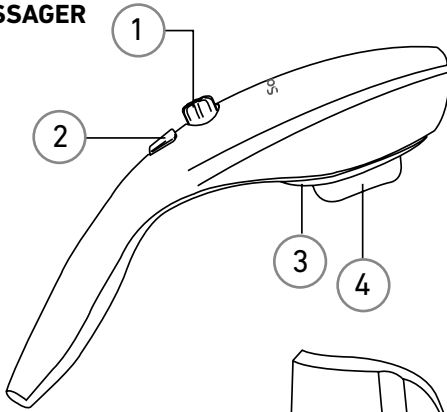1. Power/Speed knob with indicator
2. Infrared button
3. Infrared light
4. Vibration massager
5. Power port


CHARGING DOCK

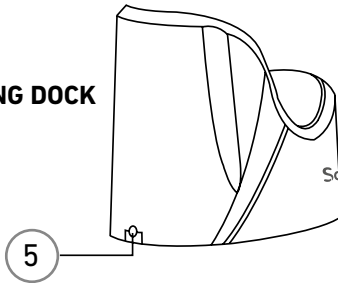
INSTRUCTIONS

POWER ON / OFF

Turn the Power / Speed knob clockwise / counterclockwise with a click sound to turn on / off the massager, the indicator will turn on / off accordingly. The default massage cycle is 15 minutes after which the massager will turn off automatically.

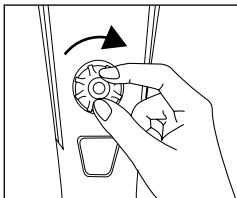
Power on

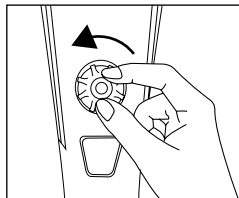
Power off

SPEED ADJUSTMENT

When the massager is turned on, turn the Power / Speed knob clockwise / counterclockwise to increase / decrease the messaging speed.


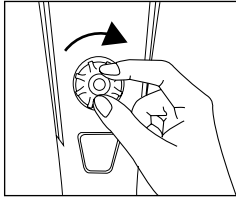
High speed

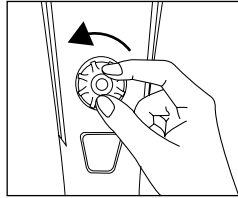Low speed

INFRARED MESSAGE

When the massager is turned on, press the Infrared Button to turn on / off infrared heat, the infrared light will turn on / off accordingly. The infrared message function can also be turned off with the power.

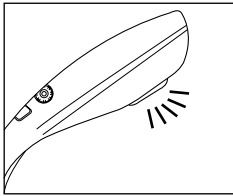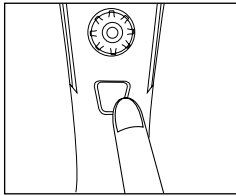CHARGING THE MASSAGER

1. Place the massager into the charging dock.
2. Plug the power cable into the power port.
3. Connect the adapter to a power source.
4. The power indicator lights up blue when charging is started.
5. The power indicator turns off automatically when fully charged after 3.5 hours.

NOTE:

- When the battery is low and the massager is turned on, the power indicator will flash.
- Do not turn on the massager when it's being charged in the dock.

NOTE


For first time users, please try the product for 10 to 15 minutes first, then extend message time if no discomfort occurs. We do not recommend using the product for over 30 minutes per session.

UNSUITABLE USERS

The following users must consult their doctor before use:

1. Under medical treatment or when not feeling well
2. Suffering from hypertension, heart diseases, brain diseases or thrombosis
3. During pregnancy or menstruation
4. Suffering from osteoporosis
5. Suffering from acute diseases
6. Having a pacemaker or other implanted devices

MASSAGE AREAS

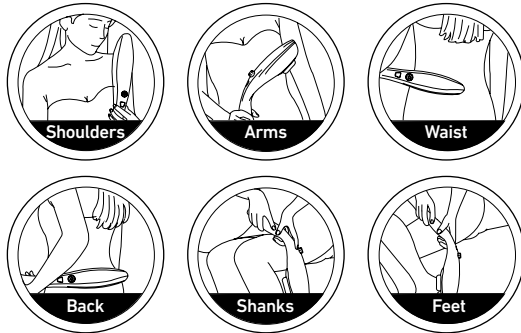TECHNICAL DATA

Input	100–240V, 50/60Hz
Output	9VDC, 1.0A, 9W
Message cycle	15 minutes
Operating time (at the highest setting)	50 minutes
Charging time	Approx. 3.5 hours

RIGHT OF COMPLAINT

By law the product must be returned to the place of purchase along with the original receipt in the event of a complaint. The user is liable for any damage to the product that is caused by not using the product for its intended purpose or by not following this user manual correctly. The right of complaint will not apply in these instances.

INFORMATION ABOUT HAZARDOUS WASTE

Used products should be handled separately from household waste. They should be taken for recycling in accordance with local rules for waste management.

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